
ordering information

Developed portfolio offering unique compatibilities such as tandem use of Catchmini & Catchmaxi inside the same 6F guiding catheter (Fargomax6f).


Catch+ mini



Catch+

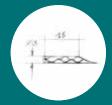

Catch+ maxi

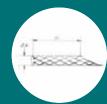

Reference	Stent Ø		Characteristics	Unconstrained stent		Compatible delivery
	Ø Nominal (mm)	Length (mm)		Ø (mm)	Length (mm)	microcatheter
CATCH MINI 3x15	3	15	For vessels up to 3 mm	4	15	VASCO+10 ID min .017''
CATCH	4	20	For vessels up to 4 mm	5	20	VASCO+18 ID min .021''
CATCH MAXI 6x30	6	30	For vessels up to 6 mm	6,5	30	VASCO+21 ID min .024''

Reference: 1. superior technology compared to the 1st CATCH version. Jeffrey L. Saver et al., Stent Retriever Thrombectomy after Intravenous t-PA vs. t-PA alone in Stroke, The NEJM.org, 2015

CATCH+ is designed for use in the flow restoration of patients with ischemic stroke due to large intracranial vessel occlusion. It is indicated to restore blood flow in the neurovasculature of patients who are ineligible for intravenous tissue plasminogen activator (IV tPA), who fail IV tPA therapy or as a supplement treatment of initiated IV tPA therapy. VASCO+ reinforced micro-catheters are intended: for injection of diagnostic or therapeutic products; to position pushable coils "SPIRALES" or detachable coils especially the ones of MDS « mechanical detachment system »; for the use of the selfexpanding stent LEO+ or SILK+. FARGOMAX are guiding catheters intended to facilitate the introduction of micro-catheters for therapeutic and diagnostic use. CATCH+, VASCO+ and FARGOMAX are class III CEO297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BAIT Extrusion. Carefully read the instructions for use before use. Not reimbursed. First CE marking: 2012 (CATCH+), 2004 (VASCO+) & 2009 (FARGOMAX).

The content of this document, in particular data, information, trademarks and logos is BALT S.A.S and affiliates' sole property ©2018 BALT S.A.S and affiliates all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT S.A.S and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes.

DC013GB (08/2018)

Bali

10, rue de la Croix Vigneron, 95160 Montmorency, Franc

Tél.: +33 (0)1 39 89 46 41 Fax: +33 (0)1 34 17 03 46

www.balt.fr



catch+

revascularization device




catch+

Designed for use in the flow restoration of patients with ischemic stroke due to large intracranial vessel occlusion.

superior technology

that significantly benefits the patient

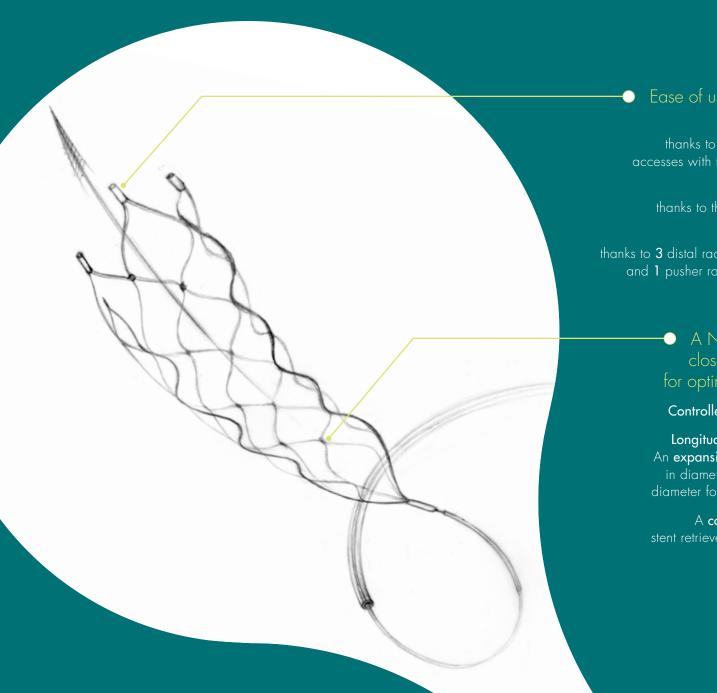

Easy navigation

thanks to proximal & distal accesses with the Vasco+ range

Easy control

thanks to the pusher stability

Easy positioning

thanks to 3 distal radiopaque markers and 1 pusher radiopaque marker

A Nitinol laser-cut

Controlled foreshortening

Longitudinal slit to allow: An expansion up to 6,5 mm in diameter (unconstrained diameter for the Maxi range)

A compression of the stent retriever in small vessels



